

Dental Suction Unit

VC10, VC30, VC30+, VC60+, VC80+

Auto-control technology: Highly efficient suction separation



VC30+

VC60+

VC80+



VC10

VC30



Efficient separation of wastewater



Two-stage air/water separation process, extra care for safety



Auto start-stop system: Reactive to the start & stop of dental units




Detecting dental unit, automatic start and stop



Double motors for vacuum pump and drainage



Low noise level

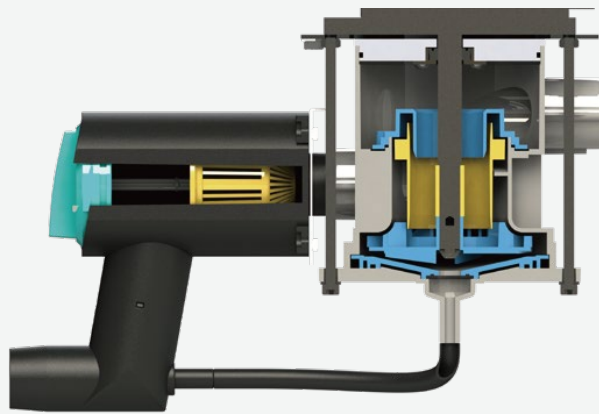

Efficient separation of wastewater

Product Details

Mode	VC10	VC30	VC30+	VC60+	VC80+
Applicable workplace	1-2	3-4	3-4	5-8	8-10
Input power	220V(1~)/50-60Hz	220V(1~)/50-60Hz	220V(1~)/50-60Hz	220V(1~)/50-60Hz	220V(1~)/50-60Hz
Rated power	500W	800W	1100W	1850W	2350W
Maximum current	4.8A	6A	7A	8A	10A
Negative pressure	-13	-17	-19	-21	-22
Max.Gas flow	1000L/min	1300L/min	1300L/min	2200L/min	3600L/min
Frequency conversion	NO	NO	NO	YES	YES
Dual-motor	NO	NO	YES	YES	YES
Noise	56dB(A)	56dB(A)	58dB(A)	62dB(A)	64dB(A)
Weight	17kg	23kg	38kg	45kg	56kg
Dimension	300*440*580mm	350*500*620mm	453-370*684mm	463*415*715mm	475-468*780mm